



LOW RADIATION DOSE

- Patients are exposed to a **low dose** pre-op **CT scan**, resulting in a lower radiation exposure than a single full spine x-ray
- Pre-operative planning potentially **eliminates the need of intra-operative checks**, with a dramatic reduction of irradiation **-33%** compared to the free-hand technique^[11,13]
- The **cumulative dose is potentially reduced** with respect to navigation-assisted technique

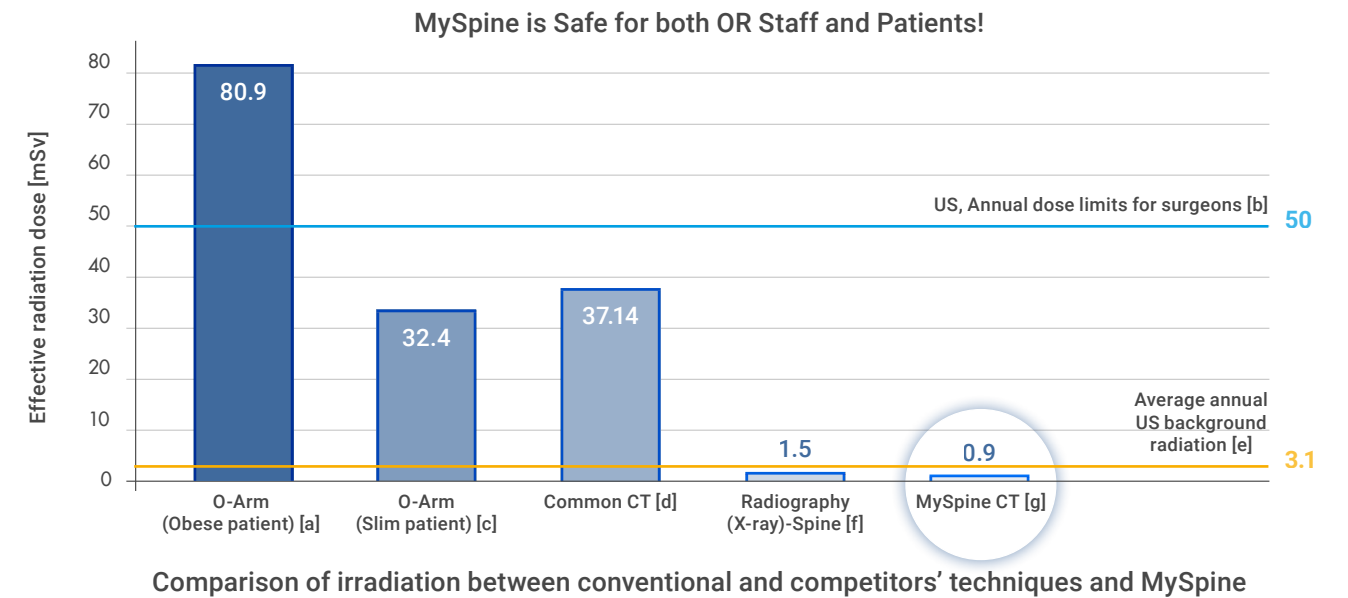

[a] Lange et.al. Estimating the effective radiation dose imparted to patients by intraoperative cone-beam computed tomography in toracolumbar spinal surgery, Spine 2013 [b] US Nuclear Regulatory Commissions (USNRC) [c] Lange et.al. Estimating the effective radiation dose imparted to patients by intraoperative cone-beam computed tomography in toracolumbar spinal surgery, Spine 2013 [d] Biswas et.al. Radiation Exposure from Musculoskeletal Computerized Tomographic Scans, JBJS Am. 2009 [e] Health Physics Society Specialists in Radiation Safety, Lawrence Berkeley National Laboratory, Fact Sheet 2010 [f] Radiation Dose in X-Ray and CT Exams; 2013 Radiological Society of North America, Inc [g] MySpine, Charité University Hospital, Berlin, Germany

MYSPINE CASE MANAGEMENT

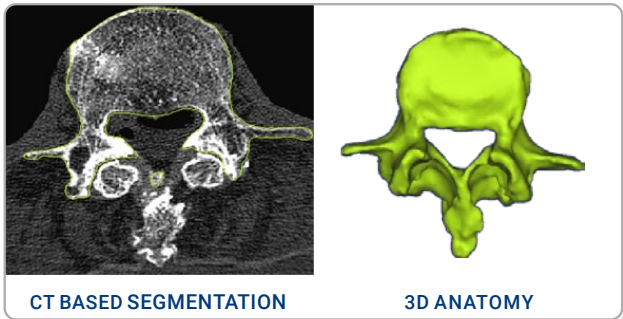

1. IMAGE ACQUISITION

Low Dose CT scan to deliver 3D reconstruction of **each patient's** vertebral anatomy




3. 3D PRINTING MYSPINE MC

3D patient-matched guides are sent to the hospital



2. 3D PRE-OP PLAN MANAGEMENT

The surgeon defines optimal implant parameters: screw **diameter**, **length** and **trajectory**

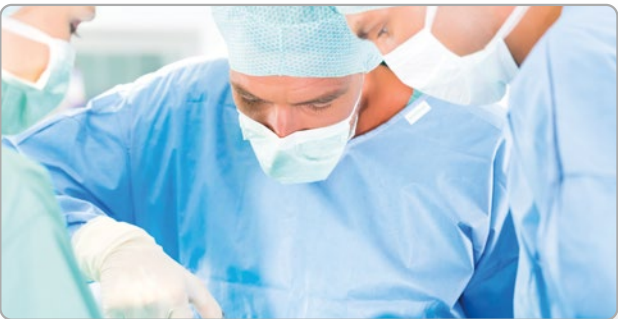

4. MYSPINE MC MIS SURGERY

Surgery with dedicated MySpine MC system

REFERENCES

- [1] Matsukawa K. et al., Incidence and Risk Factors of Adjacent Cranial Facet Joint Violation Following Pedicle Screw Insertion Using Cortical Bone Trajectory Technique, Spine, 2016
[2] Sakaura H. et al., Posterior lumbar interbody fusion with cortical bone trajectory screw fixation versus posterior lumbar interbody fusion using traditional pedicle screw fixation for degenerative lumbar spondylolisthesis: a comparative study, JNS, 2016
[3] Khanna N. et al., Spine (Phila Pa 1976). 2016 Apr;41 Suppl 8:S90-6. doi: 10.1097/BRS.0000000000001475. Medialized, Muscle-Splitting Approach for Posterior Lumbar Interbody Fusion: Technique and Multicenter Perioperative Results
[4] Gautschi O. et al., Maximal access surgery for posterior lumbar inter body fusion (PLIF) with divergent, cortical bone trajectory (CBT) pedicle-screws: a good option for minimize spine access and maximize the field for nerve decompression, Journal of neurosurgical sciences, 2015
[5] Matsukawa -2nd MORE Japan MySpine cortical Bone Trajectory 2017. <https://mysurgeon.medacta.com/uploads/presentation/attachments/d33a45ed-c550-438b-96b8-5e3fb1696725.mp4>
[6] Matsukawa - Biomechanics of CBT (internal file)
[7] Regev G et al., Nerve injury to the posterior rami medial branch during the insertion of pedicle screws: comparison of mini-open versus percutaneous pedicle screw insertion techniques. Spine. 2009;34(11):239-42
[8] Lamartina C. et al., Pedicle screw placement accuracy in thoracic and lumbar spinal surgery with a patient-matched targeting guide: a cadaveric study, European Spine Journal, 2015
[9] Santoni B.G. et al., Cortical bone trajectory for lumbar pedicle screws, The Spine Journal, 2009
[10] Mori K. et al., Short-Term Clinical Result of Cortical Bone Trajectory Technique for the Treatment of Degenerative Lumbar Spondylolisthesis with More than 1-Year Follow-Up, Asian Spine Journal, 2016
[11] Farshad M. et al., Accuracy of patient-specific template-guided vs. free-hand fluoroscopically controlled pedicle screw placement in the thoracic and lumbar spine: a randomized cadaveric study, European Spine Journal, 2017
[12] Chin K.R., Clinical Outcomes With Midline Cortical Bone Trajectory Pedicle Screws Versus Traditional Pedicle Screws in Moving Lumbar Fusions From Hospitals to Outpatient Surgery Centers, Clinical Spine Surgery, 2017
[13] Petrone S. et al., Cortical bone trajectory technique's outcomes and procedures for posterior lumbar fusion: A retrospective study, Journal of Clinical Neuroscience, 2020
[14] Matsukawa K. et al., Accuracy of cortical bone trajectory screw placement using patient-specific template guide system, Neurosurgical Review, 2019
[15] Matsukawa K. et al., Cortical pedicle screw trajectory technique using 3D printed patient-specific-guide, M.O.R.E. Journal, 2018
[16] Marengo N. et al., Cortical Bone Trajectory Screw Placement Accuracy with a Patient-Matched 3-Dimensional Printed Guide in Lumbar Spinal Surgery: A Clinical Study, WORLD NEUROSURGERY, 2019
[17] Marengo N. et al., Cortical Bone Trajectory Screws in Posterior Lumbar Interbody Fusion: Minimally Invasive Surgery for Maximal Muscle Sparing—A Prospective Comparative Study with the Traditional Open Technique, Clinical Study, 2018

All trademarks and registered trademarks are the property of their respective owners.
This document is intended for the US market.



Medacta International
Strada Regina, 34 - 6874 Castel San Pietro - Switzerland
Phone +41 91 696 60 60 - Fax +41 91 696 60 66
Info@medacta.ch - www.medacta.com

Medacta USA
6640 Carothers Pkwy - Franklin, TN 37067
Phone +1 866 830 1063 / +1 615 622 4715 - Fax +1 312 896 9138
info@medacta.us.com

MySpine® MC Leaflet
99.my46MC.11US
rev. 02
Last update: August 2023



MINIMALLY INVASIVE PATIENT-MATCHED SOLUTIONS



Brochure


Joint

Spine

Sports Med


ACCURATE PRE-OP PLANNING

The MySpine Web Platform allows for a **simple** and **accurate 3D pre-operative planning**. The surgeon can simulate the final screw position in the patient's medical images and **preview any potential surgical critical aspects**.

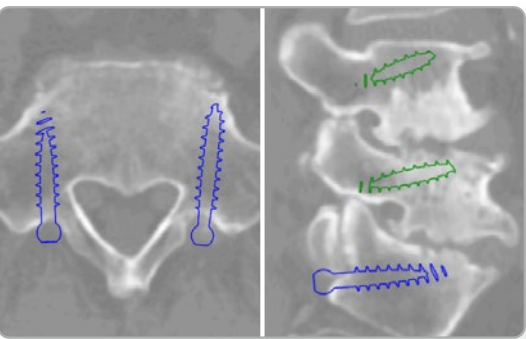

2D CT VIEW

An **effective** tool for a **personalized surgical planning**.

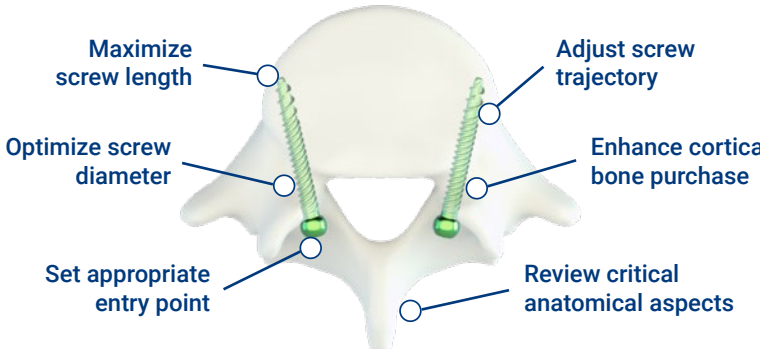

3D VIEW

- Maximize screw length
- Adjust screw trajectory
- Optimize screw diameter
- Enhance cortical bone purchase
- Set appropriate entry point
- Review critical anatomical aspects

MINIMALLY INVASIVE PATIENT-MATCHED SOLUTIONS

MySpine MC is a **3D printed** patient-matched solution in the **midline cortical** approach. Posterior lumbar fusion is performed in a **minimally invasive**, muscle sparing way, allowing for shorter operating times and a substantial reduction of both radiation exposure and costs.

- **MINIMALLY INVASIVE**
- **EXCELLENT CLINICAL OUTCOMES**
- **HIGH EFFICIENCY**
- **LOW RADIATION DOSE**

The goal of MySpine MC is to combine an **excellent fusion rate** with greater **predictability** of the clinical outcomes.





MySpine **MC** - Midline **C**ortical

EXCELLENT CLINICAL OUTCOMES

MySpine MC provides **highly precise implant positioning** which:

- allow **accurate positioning** of entry points in the the pars interarticularis with **favourable cortical bone purchase**^[4]
- may enable the use of **longer screws** and **larger diameters** than CBT free hand^[5]
- may lead to **uncompromised fusion rate**^[6]

Moreover, the pre-operative trajectory management may **reduce the risk of nerve root injury**^[7]

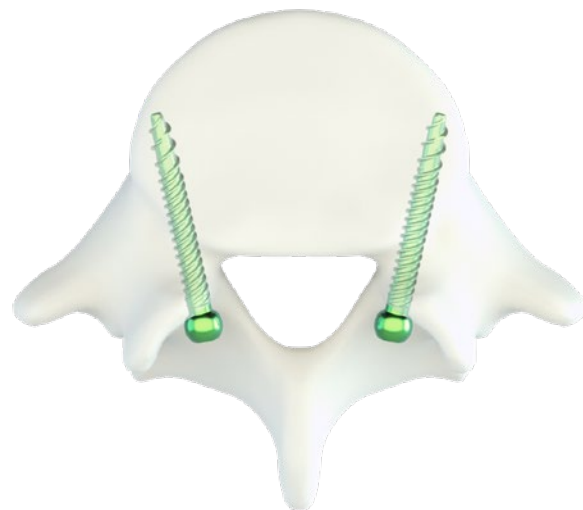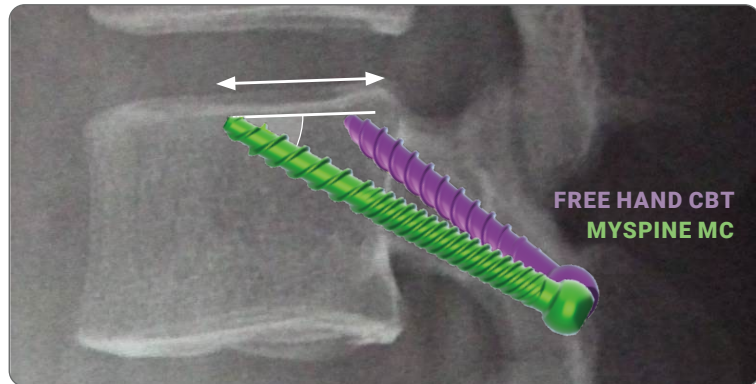

99.5%

SAFE PEDICLE SCREW
POSITIONING^[14]

-69%

REDUCED SCREW
LOOSENING RATE^[9]

+35%

SIGNIFICANT INCREASE
IN PULL-OUT RESISTANCE^[14]


-83%

ANTEROPOSTERIOR
SPONDYLOLISTHESIS
CORRECTION SLIP^[10]

WHY A MYSPINE MC MINIMALLY INVASIVE SURGERY?

Thanks to its muscle sparing technique, the erector spinae muscles are gently manipulated and a **small skin incision** of 4-5cm is performed.

For this reason, MySpine MC delivers a **minimally disruptive surgery**, which is fundamental to drive a **fast patient recovery**. MySpine MC has the potential to **improve the patients' quality of life** and **help support their recovery** after a spinal fusion surgery.

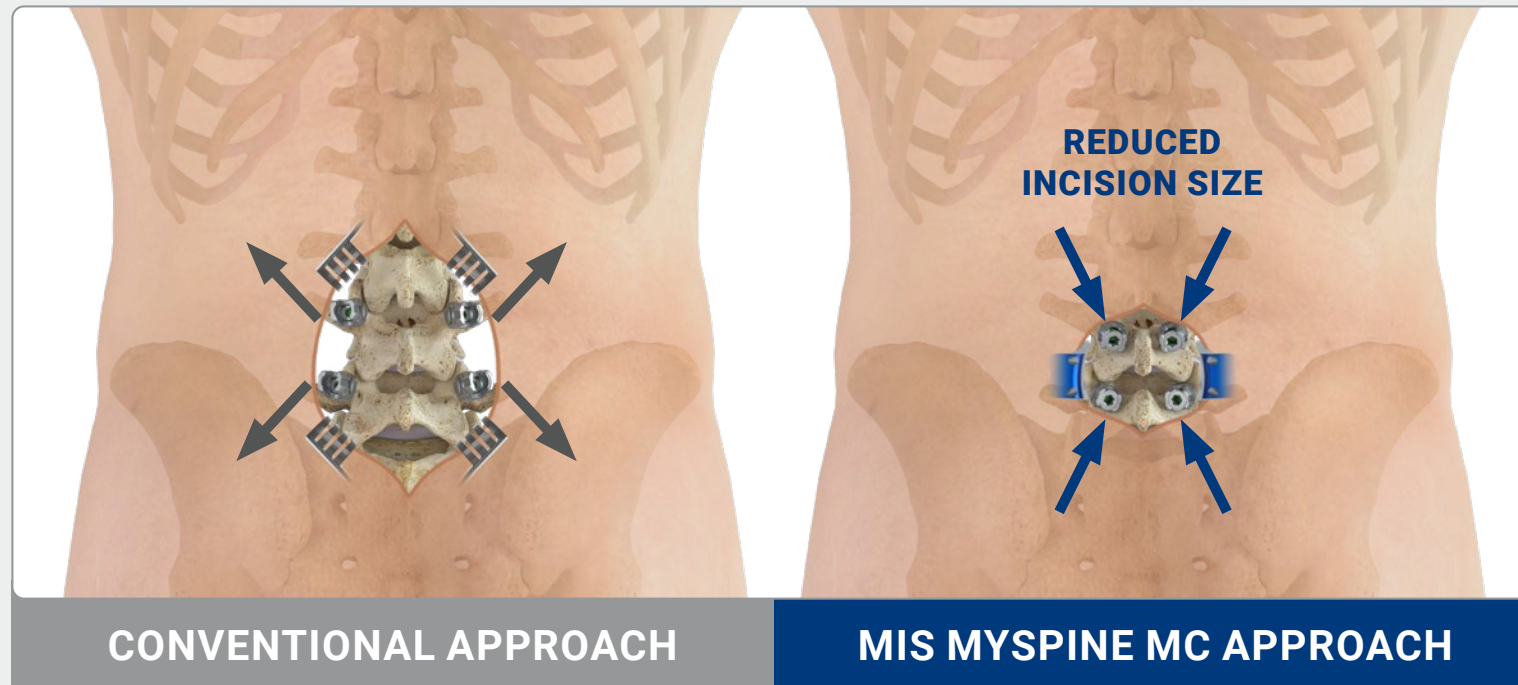

From Minimally Invasive Surgery
to Personalized Medicine
and beyond

DECREASED POST-OPERATIVE PAIN

In comparison with "conventional" open surgical techniques, the MySpine MC approach may **reduce the postoperative pain** thanks to a **less invasive technique**.^[16,17] ODI index at 12 months is **reduced by 18% more** than conventional technique, leading to a **better patient clinical score**.

-18%
POST-OP PAIN^[16,17]

SHORTER REHABILITATION

While **not violating the neuro-muscular structures**, the MySpine MC technique may **decrease the muscular atrophy** leading to a **shorter rehabilitation**.^[16,17]

"My patients can **walk autonomously** the day after the surgery." MD N. Marengo, Italy

SHORTER HOSPITAL STAY

The MySpine MC technique usually **significantly reduces the duration of the hospital stay by 37%**.^[13]

"MySpine MC is a **Minimally Invasive** technique proven to be successful in **Outpatient Setting**." MD I. LaMotta, USA

-37%
HOSPITAL STAY^[13]

SMALL SKIN INCISION

With MySpine MC, the skin incision is often shorter than with "conventional" open surgery and therefore **scar tissue is reduced**, guaranteeing an **easier soft tissue handling** and a **more "cosmetic" procedure**.^[16,17]

FASTER RETURN TO DAILY ACTIVITIES

The MySpine MC 3D Printed Patient-Specific Solution may provide **better biomechanical performance**, allowing for an **improved long-term outcome**.^[15,16,17]

"At 6-month follow-up, our patients show **important clinical improvements**, without new neurologic deficits or radiologic pathologic findings." MD K. Matsukawa, Japan

LESS BLOOD LOSS

Preservation of muscles and vessels may **reduce blood loss**, -16% compared to conventional open access surgery, for **more conservative treatments**.^[16,17]

-16%
BLOOD LOSS^[16,17]

REDUCED COMPLICATIONS


The MySpine MC technique **significantly reduces the incidence of complications**, when compared to free-hand techniques, because of the **highly accurate implant positioning**.^[13,14]

"In our specific setting, the same surgical team **reduced complications** from 16% using the free-hand technique to **0% with MySpine MC**." MD S. Petrone, MD N. Marengo et al., Italy

HIGH EFFICIENCY

ECONOMIC EFFICIENCY

- **No expensive** capital investment is required
- **No recurring service cost**
- **Rapid Learning Curve** for effective accuracy
- **Outpatient Surgery**: Hospitals and Ambulatory Surgical Centers can potentially capitalize on resources and increase volumes as patients return home the same day of the surgery^[12]

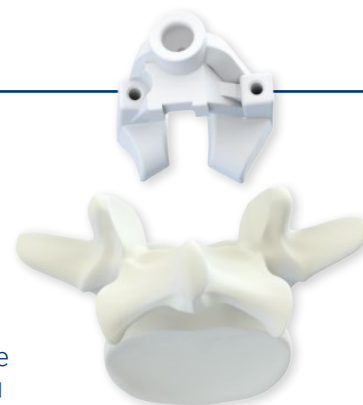

TIME EFFICIENCY

- **No peri-operative image acquisition**, thanks to the accurate pre-op planning^[11]
- Compared to free-hand CBT, the MIS MySpine MC technique leads to a significant **34% reduction of procedural time**.^[13]

-34%
PROCEDURAL TIME^[13]