

MySpine® CERVICAL

PATIENT MATCHED TECHNOLOGY
IN SPINE SURGERY

UNIQUE ANATOMIES, PATIENT-MATCHED SOLUTIONS



Brochure

Joint

Spine

Sports Med

MySpine Cervical is a patient-matched technology to be used as **anatomical perforating guides**, specific to a single patient's anatomy, for the intraoperative preparation of the pilot hole for **screw positioning in cervical spine**. The MySpine platform uses the patient's spinal CT scan to create a 3D preoperative plan specific to that patient.

The MySpine Cervical guides are indicated for levels from C2 to C7.



PEDICLE TRAJECTORY


- Pedicle screws have superior biomechanical properties compared to lateral mass screws, with a higher pullout strength^[1]
- Pedicle screws are ideal for strong fixation in cases of deformity, trauma or tumors^[1]
- Significantly longer screws compared to lateral mass trajectory^[2]



GUIDE DESIGN

- The MySpine Cervical guides are available in BILATERAL and MONOLATERAL configuration
- The spinous process blocks are available in open, open cranial, open caudal or closed option
- The monolateral guides help reduce muscle retraction while ensuring the same accuracy and safety as the bilateral guides^[3]

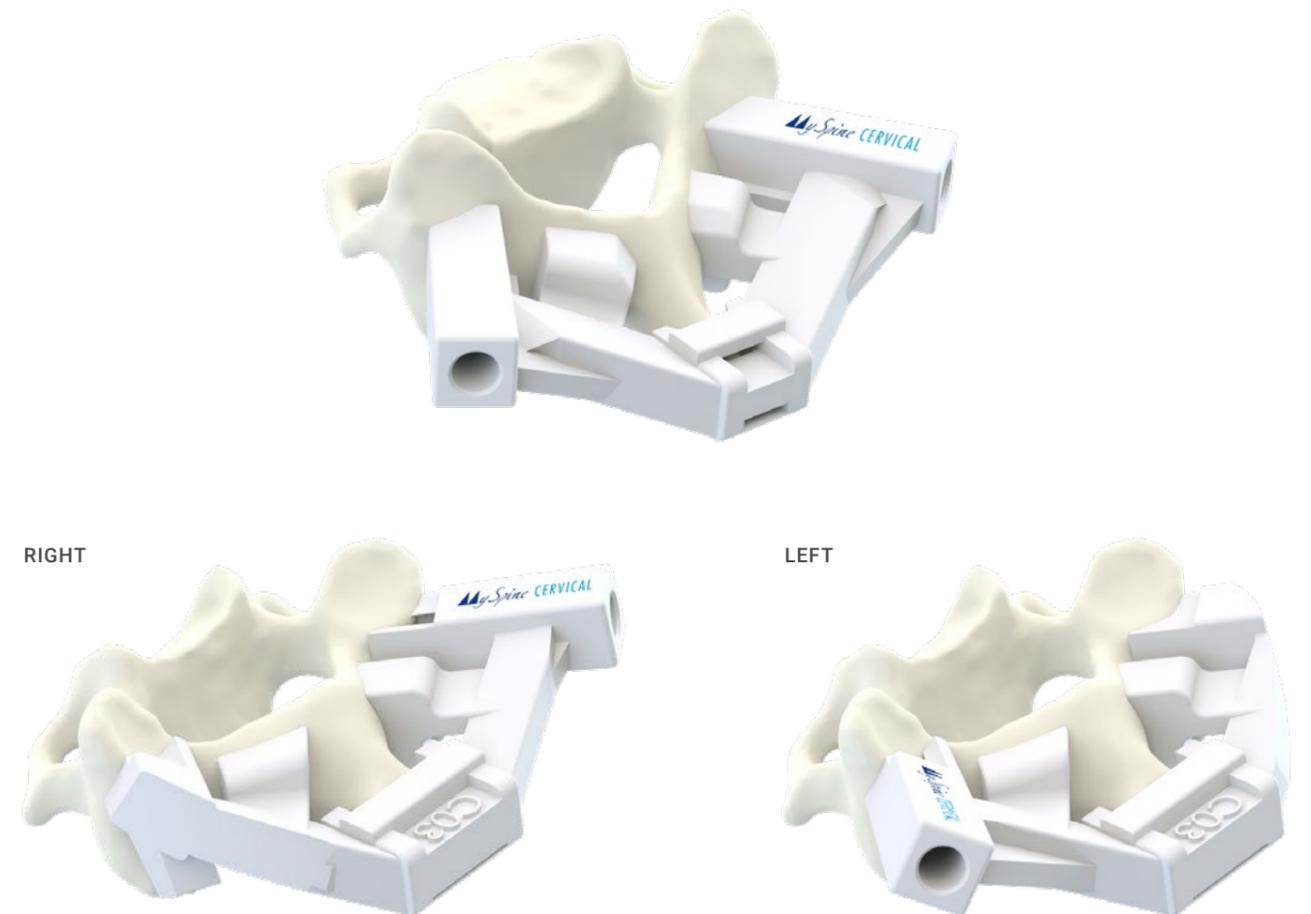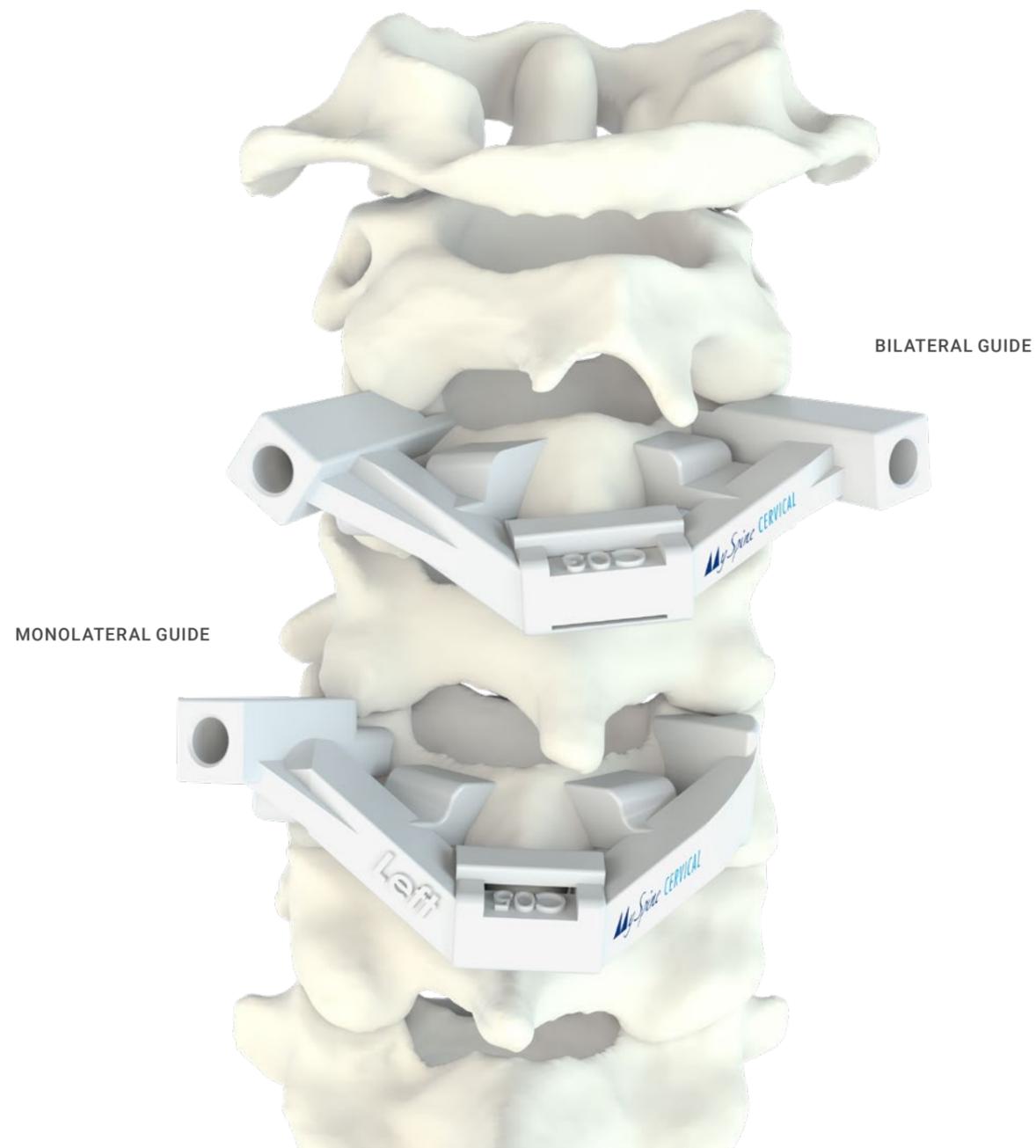

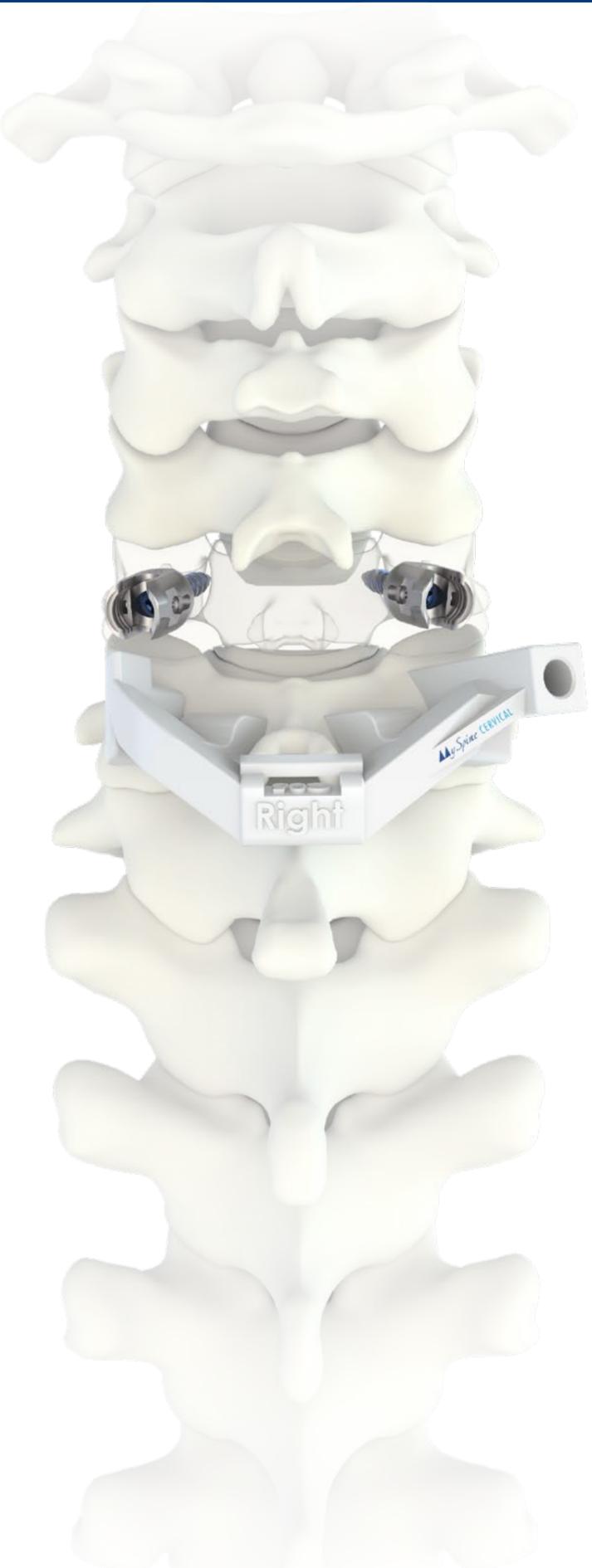

VALUE PROPOSITION

98.7% screw accuracy

MySpine Cervical allows for precise screw insertion in the pedicle, while reducing the incidence of complications.^[4]

100% safety

The precision and safety of MySpine Cervical make it a valuable navigation tool for pedicle screw insertion in the cervical spine.^[5]

49.8% reduction in O.R. Time

The average operative time with MySpine Cervical is 168 minutes, compared to a range of 110 - 564 minutes with the free-hand technique.^[5]

41.4% reduction in blood loss

The median blood loss using MySpine Cervical is 300 mL, compared to a range of 150 - 1300 mL with the free-hand technique.^[4]

Low-dose exposure

The median intraoperative fluoroscopic dose has been reduced to 321.2 mGycm² vs. the standard range of 102.4 to 825.0 mGycm².^[5]

Treatment of tumor cases


MySpine Cervical could increase the safety and the accuracy of this challenging screw placement, even in oncological patients.^[6]

CASE MANAGEMENT


Accurate 3D preoperative planning based on a **low-dose CT scan** delivers patient-matched guides, requiring zero capital investment.

IMAGE ACQUISITION

A specific **low-dose CT protocol** ensures safe image acquisition, reducing the amount of irradiation absorbed by the patient. Preoperative planning potentially nullifies the need for intraoperative checks, with a dramatic reduction of irradiation. The cumulative dose is potentially reduced compared to the navigation-assisted technique.

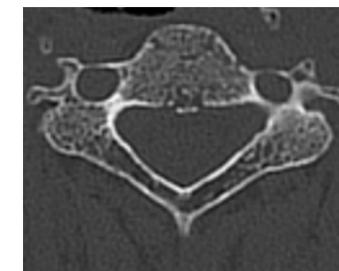


CT BASED SEGMENTATION

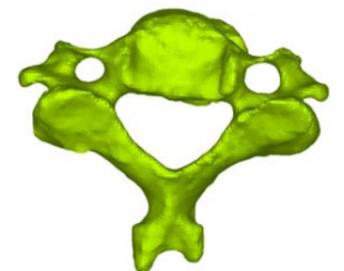

3D ANATOMY

MYSPIKE WEBPLANNER

The MySpine Cervical cases are managed entirely online, with no need to install additional software. The case database is available to the surgeon anytime and anywhere. The information on the website is always kept up-to-date.

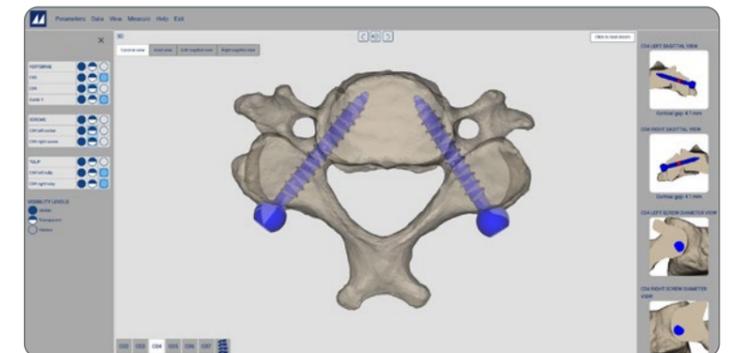

3D PRINTING

The MySpine guide is designed to match the vertebral anatomy of the patient and provide maximum stability on the vertebra, as well as correct placement of the screws. A plastic 3D model, which accurately replicates the patient's vertebra, is provided to simulate the precise placement of the MySpine guide in the operating room.



PROCTORED SURGERY

Prior to the surgery, surgeons will receive the MySpine Cervical instruments and the plastic replica of the patient vertebrae. Medacta's mission is to provide surgeons support through the M.O.R.E. Education Program. Upon request, the surgeon can receive assistance from an experienced Medacta Proctoring Surgeon to support them during their first surgeries in the local hospital or surgery center. With Medacta, the surgeon is never alone.




SYNERGY

The M.U.S.T. MINI posterior cervical screw system is a **comprehensive solution** for fixation of the occipito-cervico-thoracic spine. The variety of screws, hooks, rods and connectors allows the surgeon to tailor the construct to the specific patient anatomy and pathology to be treated.


High SAFETY and ACCURACY in cervical screw placement^[4,5] with

M.U.S.T. MINI POSTERIOR CERVICAL SCREW SYSTEM

SCREWS

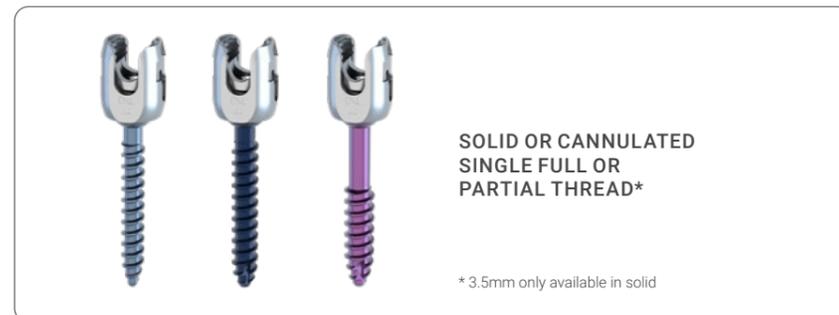

OCCIPITAL PLATE

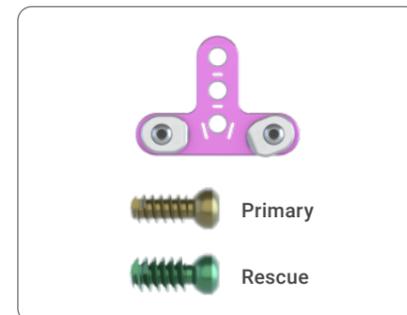

CONNECTORS & HOOKS



RODS



MYSOLUTIONS PERSONALIZED ECOSYSTEM



A network of advanced digital solutions designed to improve patient outcomes and healthcare efficiency

Medacta's MySolutions Personalized Ecosystem is designed around the patient's needs and expectations, in collaboration with an international network of expert surgeons, with the aim of delivering value throughout the entire patient journey. Surgeons' advanced 3D planning is at the core of our platform, complemented by high-precision execution tools such as patient-matched surgical guides, as well as an augmented-reality-based surgical platform and verification software. Medacta has developed a pathway to improve patients' surgical experience and support them throughout the continuum of care.

Why choose Medacta personalized enabling technologies?

PERSONALIZED 3D PLANNING

Leverage Medacta's surgeon experience with our personalized **3D planning tool**. The preoperative planning software offers solutions to accurately plan even the most challenging case.

PRECISE EXECUTION

Provide **complete and precise pre- and intraoperative guidance** for screw placement, allowing surgeons a personalized approach based on each **patient's unique anatomy**.

STREAMLINED WORKFLOW

The surgeon's procedural workflow can be seamlessly integrated with Medacta's personalized enabling technology, potentially **saving O.R. space, setup time, and increasing operating efficiency**.

FULLY SUPPORTED ADOPTION

The **learning curve** can be limited to a few cases with the support of tailored **high-level educational pathways**. With the **M.O.R.E. Institute**, the surgeon is never alone when discovering new technologies.

SUSTAINABLE SOLUTION

Limited capital investment required.





**REDEFINING BETTER
IN ORTHOPAEDICS
AND SPINE SURGERY**

MEDACTA.COM



Medacta International SA
Strada Regina 34 - 6874 Castel San Pietro - Switzerland
Phone +41 91 696 60 60 - Fax +41 91 696 60 66
info@medacta.ch

Find your local dealer at: [medacta.com/locations](https://www.medacta.com/locations)

All trademarks are property of their respective owners and are registered at least in Switzerland.
This document is not intended for the US market.

For all indications, contraindications, warnings, precautions, and potential adverse effects of the product presented, see the instructions for use or contact your local representative. For additional product information, visit www.medacta.com.

MySpine® Cervical
Brochure

ref: 99.MY46C.11
rev. 00

Last update: March 2025

References: [1] Mohamed E. et al., Lateral Mass Fixation in Subaxial Cervical Spine: Anatomic Review. *Global Spine Journal* 2021. [2] Marengo N et al. Cervical pedicle screw placement with patient-specific 3D-printed guides: accuracy and safety in a clinical experience. *Eur Spine J.* 2025. [3] Ito Z. et al. Pedicle Screws Can be 4 Times Stronger Than Lateral Mass Screws for Insertion in the Midcervical Spine. *Journal of Spinal Disorders and Techniques.* 2014. [4] Oda I, Kawaguchi Y, et al. Accuracy of pedicle screw placement using patient-specific template guide system. *J Orthopaedic Science.* 2021 Jan. [5] Farshad M. et al. Template guided cervical pedicle screw instrumentation. *North American Spine Society Journal.* 2022. [6] Marengo N. et al., 3D-printed guides for cervical pedicle screw placement in primary spine tumor: Case report and technical description. *Frontiers in Surgery.* 2022.